



UNIVERSITÉ GRENOBLE ALPES  
HÁSKÓLA ÍSLANDS  
UNIVERSITETET I BERGEN

Baptiste Bordet (Baptiste.Bordet1@protonmail.com)

---

# Documentation of Automatic Liquid Water Isotope Calibration-tool (ALWIC-tool)

---



LAST UPDATE : 28/02/2023

*Project in collaboration with :*

Hans Christian Steen Larsen

Daniele Zannoni

Arny Sveinbjörnsdottir

---

# Contents

<b>1</b>	<b>Introduction to ALWIC-tool</b>	<b>1</b>
1.1	Memory effect . . . . .	1
1.1.1	Van Geldern correction method . . . . .	2
1.1.2	Gröning correction method . . . . .	2
1.2	Calibration . . . . .	3
<b>2</b>	<b>How to install ALWIC-tool ?</b>	<b>3</b>
2.1	Download from github . . . . .	3
2.2	Install ALWIC-tool with Anaconda . . . . .	3
2.3	Install ALWIC-tool with another Python distribution . . . . .	4
<b>3</b>	<b>How to make a run ?</b>	<b>4</b>
<b>4</b>	<b>Output of ALWIC-tool</b>	<b>6</b>
4.1	New windows . . . . .	6
4.2	Saved Results . . . . .	8
<b>5</b>	<b>Options of treatment</b>	<b>8</b>
5.1	Using google drive . . . . .	9
5.2	$\delta^{17}O$ mode . . . . .	9
5.3	Control samples . . . . .	9
5.4	Prefill . . . . .	9
<b>6</b>	<b>Miscellaneous</b>	<b>9</b>
6.1	Flags . . . . .	9
6.2	Set-up a google account to interact with ALWIC-tool . . . . .	10
6.3	Set-up the prefill function . . . . .	11
6.4	Standards true values . . . . .	11
6.5	Exponential parameters . . . . .	12
6.6	Changes in the code . . . . .	12
<b>7</b>	<b>Future developments</b>	<b>13</b>
7.1	Executable file . . . . .	13
7.2	Other instruments . . . . .	13
7.3	Report . . . . .	13
7.4	Uncertainty . . . . .	13
7.5	Drift correction . . . . .	13
<b>8</b>	<b>Licence</b>	<b>13</b>
<b>9</b>	<b>Acknowledgment</b>	<b>14</b>
<b>10</b>	<b>References</b>	<b>14</b>

## 1 Introduction to ALWIC-tool

ALWIC-tool is a program made to process water isotopes measurements in liquid injection mode. ALWIC-tool aims to provide an highly customisable tool while having high constraints on the user inputs to avoid as much as possible human errors. We aim this code to be suitable for a wide range of people no matter what is your knowledge on both water stable isotopes measurements or programming. We will start this documentation by a brief description of the effects occurring during measurements of water isotopes.

### 1.1 Memory effect

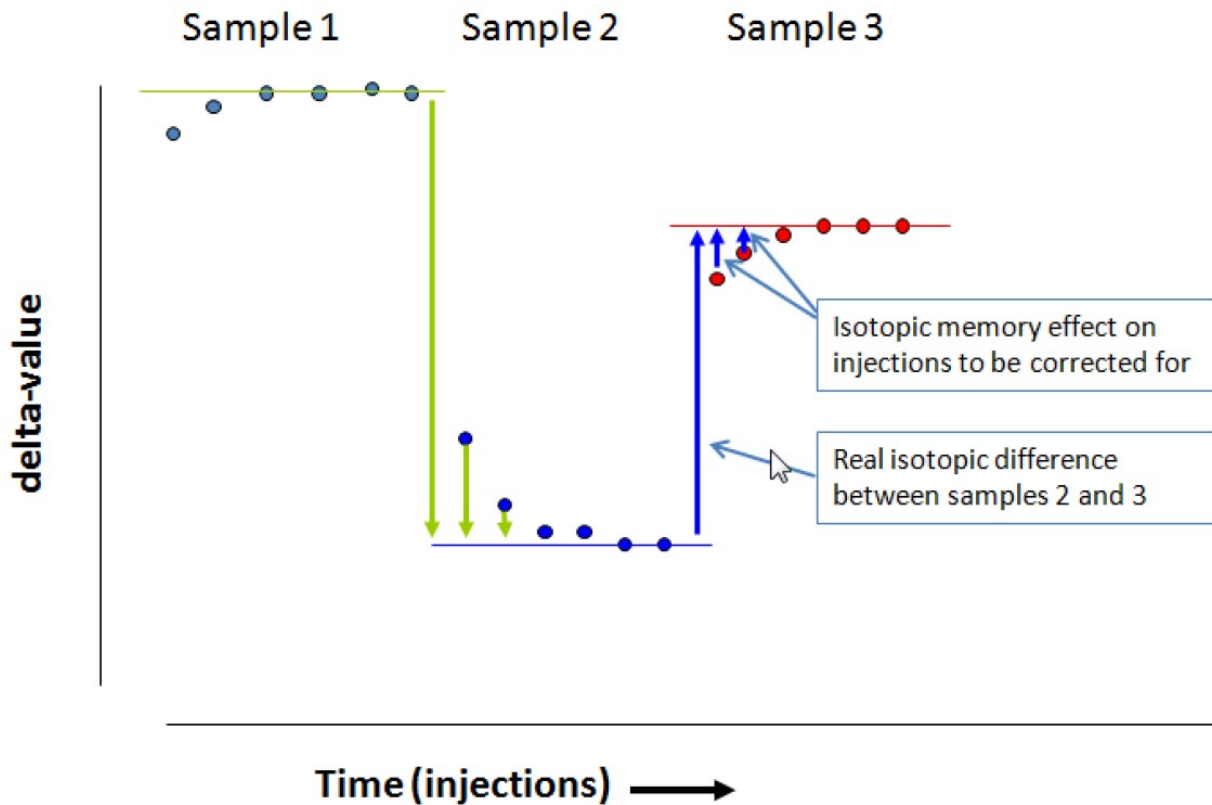


Figure 1: Illustration of the memory effect, taken from the SICalib user manual (Gröning [2018])

Since the injection system of water isotopes analyser is not perfectly hydrophobic some water will stick to the walls of the system (tubing, syringe mostly) during the injection of the sample. When injecting the next sample, some of this residual water will be blown away from the injection system and will be measured along with the new sample. This effect called memory effect produce results as shown in figure 1. The first injections of one sample are influenced by the value of the previous sample and this influence decreases with the number of injection. In ALWIC-tool there are two methods implemented to correct this effect:

- The method suggested by van Geldern and Barth [2012] who proposed to evaluate the proportion of residual water in the mixing with standards and correct on samples.

- The method suggested by Gröning [2011] who proposed to evaluate the residual water through a two memory reservoir.

### 1.1.1 Van Geldern correction method

To evaluate the proportion of residual water, a constrained non-linear optimisation is applied on standards' injections. This only means to adjust the memory coefficients (i.e. the proportion of water at each injection) to minimize the standard deviation of all the standards combined. The method used to calculate the memory coefficients is a constraint non linear optimization performed with the COBYLA algorithm (Constrained Optimization BY Linear Approximation). Details on the implementation can be found here : Scipy Community [a], Scipy Community [b].

From these coefficients the correction of all the values can be made through the formula given in van Geldern and Barth [2012] :

$$\delta_{j(mc)}^n = \delta_{j(raw)}^n + (1 - m_j) * (\delta_{j(raw)}^n - \delta_t^{n-1}) \quad (1)$$

where :

- $\delta_{j(mc)}^n$  is the memory corrected value of the  $j^{th}$  injection of the  $n^{th}$  sample,
- $m_j$  is the memory coefficient of the  $j^{th}$  injection
- $\delta_j^{n-1}$  is the true isotopic value of the previous sample.

### 1.1.2 Gröning correction method

The method is subdivided into two steps :

- First, there is the single factor calculation which calculates an approximation using the last three injection to correct the current injection.
- Then, an exponential where the parameters are set manually (Will be set to automatic if it's possible) is used to correct more accurately each injection.

Gröning's method is based on two successive correction. The first one assumes a single memory reservoir with "a constant relative memory decrease with time" (Gröning [2018]). The second one assumes two memory reservoir represented by two exponential laws combined to represent their interactions.

The first part is all done automatically. It's based on the evaluation of the memory contribution  $f$  for each pair of samples with the following formula :

$$f = \frac{a_j}{a_{j-1}} = \frac{true_j - meas_j}{true_j - true_{j-1}}$$

$j$  refer to the sample number and starts at the second sample.

To take into account of the last three samples  $f$  is reduced :

$$f = f * (1 - f - f * f)$$

Then a mean of all  $f$  is done excluding the ones with two little difference in  $a_i$  ( $a_i$  should be more than 0.4 ‰ for  $\delta^{18}O$ , 4 ‰ for  $\delta D$  and 0.4 ‰ for  $\delta^{17}O$ ) and the  $f$  negative (physically impossible).

Finally the first correction is applied with the following equation :

$$\delta_i^{first\ correction} = \delta_i^{raw} + f * (\delta_i^{raw} - \delta_{i-1}^{raw}) + f^2 * (\delta_i^{raw} - \delta_{i-2}^{raw}) + f^3 * (\delta_i^{raw} - \delta_{i-3}^{raw}) \quad (2)$$

Then to start the second correction, a percentage is calculated based on the following formula :

$$corr1 = \frac{\delta_i^{first\ correction} - true_j}{true_{j-1} - true_j} * 100$$

This percentage is then adjusted with the exponential parameters, alpha, beta and balance. These parameters are specific to each instruments.

$$corr = corr1 * (balance * e^{(-\alpha*(i-1))} + (1 - balance) * e^{(-\beta*(i-1))})$$

Finally this correction term is applied to the data :

$$\delta_{mc} = \delta_i^{first\ correction} - corr * \left( \frac{true_{j-1}}{true_j} * 100 \right) \quad (3)$$

## 1.2 Calibration

The calibration is done by measuring standards with the instruments and comparing their values to the values measured by references

The values are then corrected again with the standards. The memory corrected standards (all except the first one which is used as a rinsing one) can be used to create a calibration curve. The calibration curve is evaluated with the scipy linear regression function (Scipy Community [c]).

## 2 How to install ALWIC-tool ?

### 2.1 Download from github

In order to download ALWIC you can click [here](#). This will download a zip file with all the modules. If you are familiar with github you can clone the project with the following command :

```
git clone (repo_link) <local directory>
```

### 2.2 Install ALWIC-tool with Anaconda

The first recommendation is to create a new environment to avoid any conflicts of the installed packages with your previous configuration. Open an Anaconda prompt. Use the following command to go to the unzipped folder of ALWIC-tool :

```
cd path_to_the_folder
```

To create the new environment use the command below :

```
conda env create -f requirements.yml
```

Note : This command might take a few minutes to run since all the packages required for AIWIC-tool will be downloaded and installed.

This will create an environment called ALWIC\_env which can be activated with the command :

```
conda activate ALWIC_env
```

The environment has to be activated each time you want to run ALWIC-tool. It can also be done in Anaconda Navigator by choosing it instead of base on top left of the window. After this you can open the file called ALWIC\_main.py in Spyder and run it.

## 2.3 Install ALWIC-tool with another Python distribution

You will need to install the packages included in the requirements.txt file. This file is included in the file you downloaded. The method to install it is yours, depending on your IDE but the command to install all the required packages is :

```
pip install -r requirements.txt
```

Then you can run the ALWIC\_main.py file and start processing.

## 3 How to make a run ?

The screenshot shows the ALWIC-tool interface with the following components and annotations:

- 1**: Help button (circled in red).
- 2**: Configuration section with two dropdowns: "Protocol to follow : Grönning d17O mode" and "Choose the type of directory : Local directory".
- 3**: "Welcome to ALWIC-tool !" message and a filename input field containing "HKDS2156\_IsoWater\_20230106\_162252" with a "Prefill" button.
- 4**: "Number of standard injected 4" with an "Update" button.
- 5**: "Spy samples ?" section with a checked checkbox and "Number of known samples 5" with an "Update" button.
- 6**: Injection parameters section with three rows of input fields:
 

Number of samples injected :	50
Number of injections per sample :	4
Number of injections removed per sample :	1

 and
 

Number of injections per standard :	10
Number on injections removed per standard :	1
- 7**: Operator and Processor information fields: "Operator ID : Operator" and "Processor ID : Processor".
- 8**: "TEST" table showing data for alpha, beta, and balance:
 

TEST	$\delta^{18}\text{O}$	$\delta\text{D}$	$\delta^{17}\text{O}$
alpha	1.04	1.04	1.04
beta	0.21	0.21	0.21
balance	0.7	0.7	0.7

A "Processing" button is located at the bottom right.

Figure 2: Screenshot of the interface of ALWIC-tool

### 3 HOW TO MAKE A RUN ?

---

Here will be described how to run ALWIC-tool. The details of the options presented here are further discussed in section 5.

After starting a run a first page appears similar to the figure 2. There are different parts on this page numbered from 1 to 7 :

- 1 : Help links, they will redirect either directly to this document, either to the github page of ALWIC-tool.
- 2 : Here you will choose which method to correct the memory effect. You also have to choose if you want to correct the  $\delta^{17}O$  as well. On the second box you need to chose from you want to extract the file (see section 6.2 for how to set-up a drive account).
- 3 : Here you have to write the filename without extension. If you chose the local directory option on the previous point this will be already filled. Once you done this you can use the prefill function to prefill the points 4, 5 and almost all of 6. To be able to do the prefill the file must be customized before the measurements (more information on this on section 5.4).
- 4 : Here you first need to enter the number of standard injected at the beginning of the sequence with a maximum of 6 (This can be increased on the code, see section 6.6 to do it). Then click on the update button. A table will appear where you will be able to modify only two things : the tick-boxes and the names of standards. The tick-boxes control which standards are used for the calibration (on means the standard is used for calibration). The names of the standards are controlling the true values of these standards. The values are set in a different file before running ALWIC-tool, see section 6.4 to know how to modify this file.
- 5 : Here is an optional part to fill. If you don't have any known samples in your run you can go to the next point. Else, you need to check the little tick-box and enter the number of known samples you have over the run with again a maximum of 6 known samples which can also be modified on the code, see section 6.6 for more information. Also if you're using a quality control standards you can include it here. Before clicking the update button make sure that the filename is filled because the update button will need your file to display the table. On this table you need to indicate the port where you placed these known samples. And in the in the same way as the last point you need to set-up the true values by changing the names. To modify the file controlling this you can look at section 6.4.
- 6 : Here you just have to enter a certain numbers of parameters of the run which are in order : The number of samples (with the known samples !!) injected on all the run, the number of injections per sample (must be constant along the run), the number of injections you want to remove to calculate the mean and standard deviation on each sample, the number of injection for each standard (must be the same for all standards) and the number of injections you want to remove to calculate the mean and standard deviation on each standard, see section 4 for more details on those calculations.
- 7 : Here is an optional part. You can indicate who ran the measurements and who is processing the data.

- 8 : This part is only shown when you choose one of the Gröning method. You will be able to select your instrument's name. This will then modify the parameters used for the exponential correction of Gröning. The values are set in a different file before running ALWIC-tool, see section 6.5 to know how to modify this file.

When all the steps have been made you can click on Processing. The memory correction process will start and show the first result page.

## 4 Output of ALWIC-tool

If you choose the van Geldern method to do the memory correction, a temporary window will appear before doing the calculation. This window allows you to remove either one standard if you have multiple injections with issues or you can also remove one injection (but not more) of all the standards' injections. This tool should be used carefully and only in the case of really outliers. It should not be used to improve the memory correction !

To start the processing click on next page. The window will be destroyed and the corrections will be made.

### 4.1 New windows

When the processing button is pressed and all the calculation made a new window appears. It looks like fig. 3. The parameters of the correction are shown on this page.

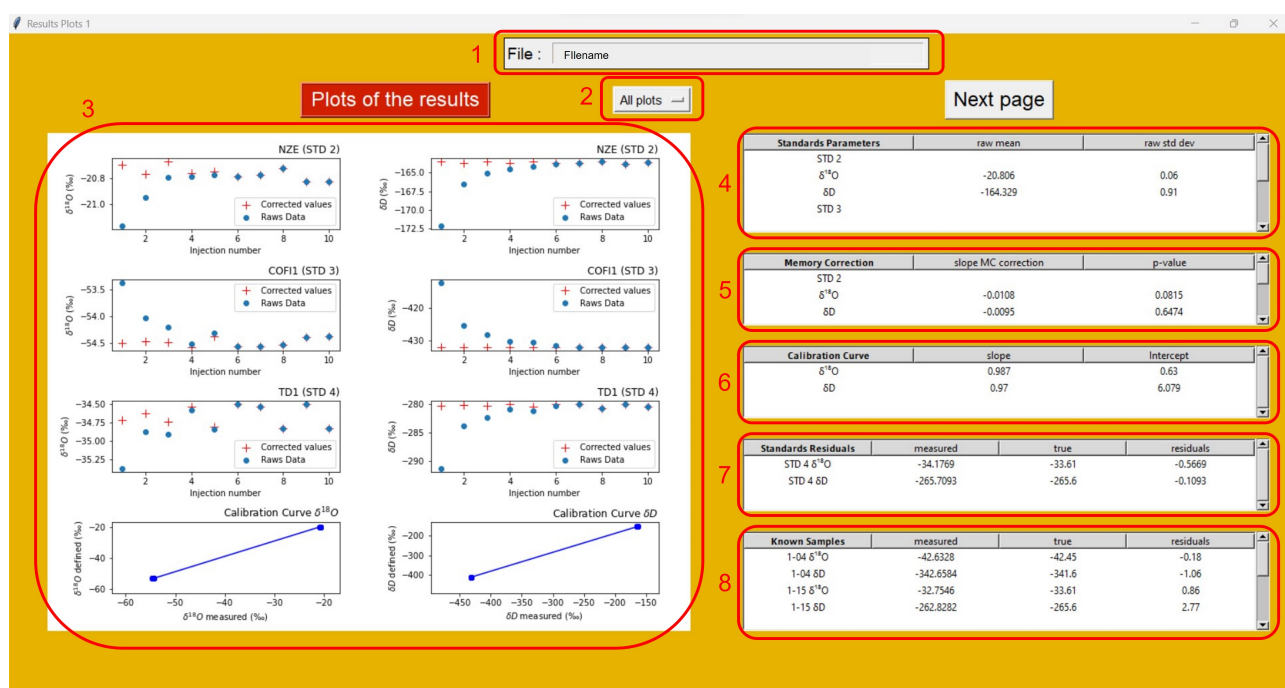


Figure 3: Screenshot of the first page of results of ALWIC-tool

As previously this page is composed of different parts from 1 to 8 :

- 1 : Here it is a reminder of the filename. It may become really useful when you start to correct several files in a row.



- 2 : When this is set to all plots the page looks like the figure 3. But if you choose an another plot the tables on the right will be overwritten to put a zoom on the plot you selected. To view again the tables, select all plots.
- 3 : You have here a summary of the correction. For each standard and isotope you will have one plot with the raw data and the memory corrected data (before calibration). These plots are auto-scaled in y so be careful before trying to compare them. At the bottom you will have the calibration curves. These plots have the same scale for the two axis.
- 4 : Here you will find the mean and standard deviation of all the standards and isotopes. The mean and standard deviation are calculated on the values you choose not to reject on the previous page. These values are calculated on the raw values.
- 5 : This table will be different according to the memory correction method you use. If you used the van Geldern method, it shows the slope of the memory corrected values and its p-value. If you used the Gröning method, it will show you the mean single factor and alpha, beta and the balance term for each isotope.

Note : Be very very careful with the p-values, it is calculated on a very limited number of values and p-value is a probability not a statistic (Wikipedia Contributors). If you are not aware enough of what are the limits do not use it !

- 6 : Here you will find the parameters of the calibration curves i.e. slope and intercept.
- 7 : This table is optional it will only show up if you don't use all the standards to show the residuals on the standard not used for calibration. The residuals are always defined by  $residuals = measured - true$ .
- 8 : Also an optional table. It will show up if you have known samples in your run. Again it will show the residuals of the known samples. The residuals are always defined by  $residuals = measured - true$ .

By clicking on next page you will not close the first page of results and will open the second page of results (figure 4). On this page the results shown are the samples values. Again there are several parts on this page from 1 to 4.

- 1 : Same as previous page, a reminder of the filename.
- 2 : This part shows several plots against the injection number.
- 3 : This table shows every mean sample (the mean being calculated by removing the number injection set at the beginning of the post-processing).
- 4 : This allows to save the data in different files (next section will describe it). You will be asked for the folder where you want to save the data. A pop-up window will appear after saving.

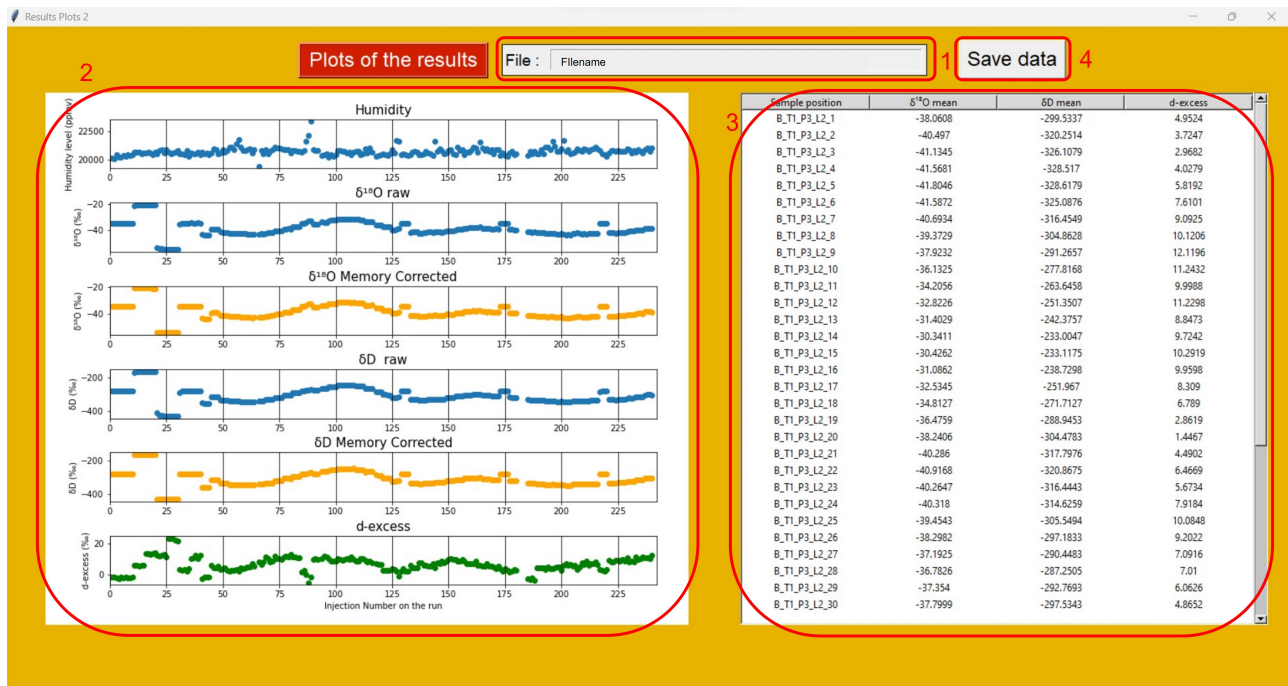


Figure 4: Screenshot of the second page of results of ALWIC-tool

## 4.2 Saved Results

There are up to 4 files after each processing.

- 1. The first file to be saved is called `your_filename_final_file`. It contains the input file with some columns added which are the memory corrected values as well as the calibrated values and some flags (more on it on section 6.1).
- 2. The second file to be saved is called `your_filename_std_parameters`. It contains a lot of information about the run like starting and ending time, parameters of the correction, the standard deviation of the humidity (detrended) along the run, names of people who ran the file and processed the data, the correction method and the mean and standard deviation of all standards.
- 3. The third file to be saved is called `your_filename_spl_results`. It contains the name of each sample as well as the mean and standard deviation of the last injections (the one you didn't choose to remove).
- 4. The last file to be saved is called `your_filename_control_samples_results` and is optional. It is only saved when you have spy samples in your run. It contains the position of each control samples with their measured values, true values and residuals.

## 5 Options of treatment

ALWIC-tool offers numerous options to process your data. They are described in the sections below.

## 5.1 Using google drive

Google drive can be used to store measured files and treated files. ALWIC-tool has a module allowing it to communicate with google drive. It uses the google API (Application Programming Interface) which allows two machines to communicate in this case ALWIC-tool and google drive. To do so the google account has to be set-up, it is described in section 6.2.

## 5.2 $\delta^{17}O$ mode

When measuring  $\delta^{17}O$  this column is corrected and a new column for 17O-excess is calculated. The 17O-excess column is calculated with the following formula according to Luz and Barkan [2005]:

$$\Delta^{17}O = \ln(\delta^{17}O + 1) - 0.528 * \ln(\delta^{18}O + 1) \quad (4)$$

## 5.3 Control samples

In order to look at the short term drift and to evaluate the performance of the run, "spy samples" can be added to the run. They are standards measured as samples (with the same number of injection and corrected in the same way as samples). It is then possible to evaluate the residuals along the run and the change of these residuals to identify potential issues.

## 5.4 Prefill

ALWIC-tool has a function to evaluate parameters directly from the file and to prefill it for the user. The filename needs to be filled and the following parameters are evaluated when it is possible : number of standards, True values of the standards, injection per standards (must be constant along one run), number of samples, injection per sample (must be constant along one run), if there are spy samples, if yes the port of each spy sample, the values of the spy samples.

Note : in order to prefill the true values of standards and spy sample the names of standards in the measurements files have to match with one standard name in the std\_values.csv file (see section 6.4)

The set-up of this function is described in section 6.3.

# 6 Miscellaneous

## 6.1 Flags

There are three flags in total. One is saved in the spl\_results file. It indicates a standard deviation of each sample higher than the threshold (By default, the thresholds are fixed to 0.06‰ for  $\delta^{18}O$ , 0.21‰ for  $\delta D$  and 0.1‰ for  $\delta^{17}O$ . See section 6.6 to know how to change these values).

The two other flags are saved in the final\_file. One is for the humidity level, if the humidity for one injection higher than 23000 ppmv or lower than 17000 a flag will be raised.

The other flag indicates if an injection has been corrected or not. There are three possible values : 0, the value hasn't been corrected either because it's the first standard either because there is NaN somewhere. 1, the value has been corrected. 2, the value has not been corrected and the final value has been replaced by the raw value. This has been set-up for the case where

the last injection of a sample is missing the next sample can't be corrected by any method. Instead of losing the value of this sample the raw values are used and the user is informed of this with a pop-up window when the second page of results is opened.

## 6.2 Set-up a google account to interact with ALWIC-tool

Setting-up the link between ALWIC-tool and a google account is a tedious task but once it's done, it can be a very powerful tool to work remotely or with co-workers.

The first thing to do is to create a google account. It is probably a good idea to use an account only dedicated to this task since the account will be set-up in a developer mode.

In google drive two folders have to be created : one for the raw files named `raws_data` (the name does not matter since the search from ALWIC-tool will look into the entire drive for the file), one for the treated files named `Treated_data` (the name can be changed on the code, see section 6.6). The files can be anywhere in the drive, the only requirements is to not have other folders with the same name.

Now, the google account will be set-up to use the google API. Before starting this be sure to be logged in with the account which will be set-up. If you have a different default account the solution is to log-out and then reconnect with the account you want to set-up in first. This will automatically change the default account.

The first thing to do is to enable the API. This can be done on the following link : <https://developers.google.com/drive/api/v3/quickstart/python>. When scrolling down you will see something similar to the figure 5. Click on Enable the API.

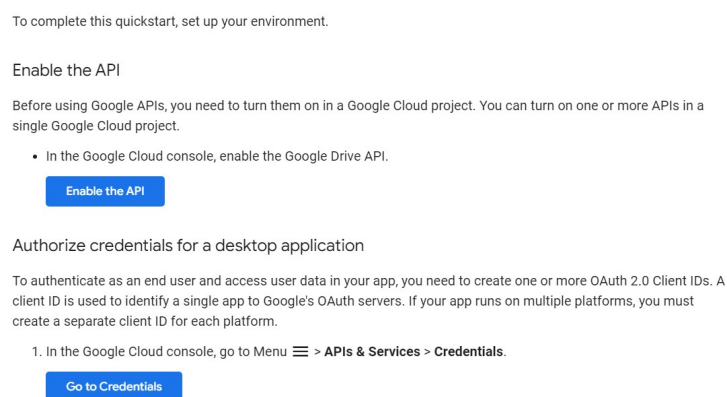


Figure 5: Box to click on to start the set-up

This will open a new window which is your google cloud console. The first thing to do is to create a new project. You chose a name and then enable the API. Then you can go back to the cloud console by clicking on Google cloud on top left of your window. Then on your left you have a panel You will have to select the credentials in API and services as shown in figure 6.

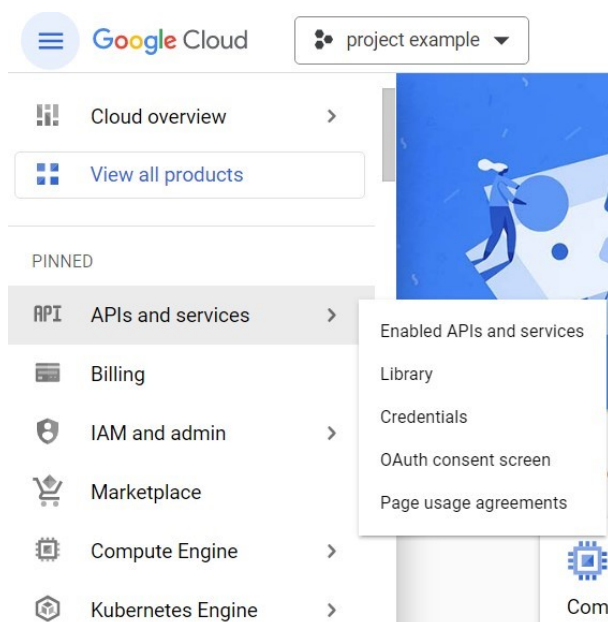


Figure 6: left panel where to find the credential control page

Then in order to create the credentials, you will have to click on create credentials and chose the API key. Then wait and close the pop up window. Then click again on create credentials and this time chose the OAuth ID. Then click on configure consent screen. Chose the user type Extern. Then, fill the App name, the user support name and the developer contact information. Skip the scope page. On the next page there the test users. You can fill up to 100 google account. Here you will fill every email which will use ALWIC-tool. This can of course be modified later. And then you can ho back to the dashboard. And click again on create crendentaials and chose again OAuth ID. Chose Desktop App as application type and then create. When the pop-up window appears download the .json file and rename it client\_code.json. You will need to move file into the app folder. If these informations are not clear enough you can refer to the video available with the code or on github. It is entitled :

### 6.3 Set-up the prefill function

In order to use the prefill function the sample descriptor file has to be set up before the measurements. If you already use it to mark your samples name you only need to add one more column with title Identifier 2 and for vial indicate if it is a standard by putting STD, a sample by putting SAMPLE or control sample by putting SPY. The code is sensitive to upper-case be careful.

If you never used the sample descriptor, there is a file example in the example\_files folder you just need to adjust it to your run.

When the Identifier 2 column is filled with STD, SAMPLE and SPY the prefill function is fully functional.

### 6.4 Standards true values

The standards values are stored in a separate file which must be closed when the code is running. The file is located in the File folder in the project folder. It is named std\_values.csv.

There are five columns in the file, the first one is the description and can be used to add some annotation on the standards. The second column is the short name, this is the name which will be displayed on screen so be careful to have meaningful names in this column. The last three columns are the values for each isotopes expressed in ‰. The file can be saved with comma or semicolon separator without distinction but the function save as should not be used. It might lead to errors in the opening of the file by python.

## 6.5 Exponential parameters

The exponential parameters are stored in a separate file which must be closed when the code is running. The file is located in the File folder in the project folder. It is named groning\_exp\_parameters.csv. The file is ordered by line. Each line is an instruments and all the parameters are ordered in different columns. The file can be saved with comma or semicolon separator without distinction but the function save as should not be used. It might lead to errors in the opening of the file by python.

The parameters can be determined with the help of SICalib (Gröning [2018]). They are very stable in time and only instrument dependent.

## 6.6 Changes in the code

- How to increase the number of standard and control sample possible in the tables ?

This requires some modifications. In the ALWIC\_main.py there are the functions called define\_table and define\_known\_sample\_table. In these function you need to follow all the comments. If you want to add standards you will also need to change a few things in the change\_table\_outliers\_page; function which is also in the ALWIC\_main below the others functions.

- How to change the defaults options when you have several modes ?

At the end of the ALWIC\_main.py file there is the definition of all the elements of GUI. If you want to change the default protocol, you will need to look at the definition of panned window 2. Changing the default folder to look for to local directory is not recommended at all because of the construction of the code.

- How to change the default values for the flag in sample result ?

In the file save\_results.py in the function save\_sample\_results\_file you can change the thresholds.

- How to change the default extension for saved file ?

At the beginning of the of the save\_results.py file, there is a variable called save\_extension which can be modified.

- How to change the Treated\_data name for google drive ?

In the file called google\_api in the function upload\_files there is a parameters query which is set by default to Treated\_data this can be changed.

Any other changes in the code can cause some unexpected side-effects. It is highly recommended to contact me via gitub (by opening an issue : (<https://docs.github.com/en/issues/tracking-your-work-with-issues/creating-an-issue>) or by creating a pull request if you want to add a feature or fix something in the code (<https://docs.github.com/en/pull-requests/collaborating-with-pull-requests/proposing-changes-to-your-work-with-pull-requests/creating-a-pull-request>)). It is also possible to send me an email but I can not guarantee I will answer you.

## 7 Future developments

### 7.1 Executable file

An executable version will be deployed soon. It will allows people to run the code without having to install any python on their machine. Stay tuned !!

### 7.2 Other instruments

For now, only Picarro's files can be treated with ALWIC-tool but in the future it will be expanded to other instruments which are doing measurements of water isotopes in liquid injection mode.

### 7.3 Report

A pdf or html document to summarize the result of each processing could be considered but for now it is too much time consuming, I will think about it only if people ask for it.

### 7.4 Uncertainty

This is definitely something that should be implemented in ALWIC-tool. However this requires a lot of time which I will not have soon. If anyone is untested in doing it, it will be a pleasure to discuss.

### 7.5 Drift correction

In agreement with Hans Christian Steen Larsen, we decided to not put a short-term drift correction because we don't think it is necessary especially with the new generation of Picarro which are more stable over time. Again if a lot of people ask for it, I will consider adding this feature.

## 8 Licence

ALWIC-tool is free software: you can redistribute it and/or modify it under the terms of the GNU General Public License as published by the Free Software Foundation, either version 3 of the License, or any later version.

ALWIC-tool is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the GNU General Public License for more details.

## 9 Acknowledgment

Thanks to Hans Christian Steen Larsen who initiated the project; Daniele Zannoni for all the support during the development of this tool; Arny Sveinbjörnsdottir for the comments made from a user perspective. Thanks to Melanie Behrens and Manfred Gröning for the help in the implementation of one of the method. Thanks to the teams at AWI, EPFL and LSCE who tested the code, your feedbacks were precious. Finally a special thanks to Pauline Gayrin and Clara Berneuil for the emotional support as well as a useful discussion.

## 10 References

Manfred Gröning. *SICalib User Manual (Stable Isotope Calibration for routine  $\delta$ -scale measurements)*. Terrestrial Environment Laboratory, International Atomic Energy Agency, Vienna, Austria, 2.16j edition, May 2018. [https://nucleus.iaea.org/sites/ReferenceMaterials/Shared%20Documents/Publications/TechnicalNotes/TELTechNote01\\_SICalibUserManual.pdf](https://nucleus.iaea.org/sites/ReferenceMaterials/Shared%20Documents/Publications/TechnicalNotes/TELTechNote01_SICalibUserManual.pdf).

Robert van Geldern and Johannes A.C. Barth. Optimization of instrument setup and post-run corrections for oxygen and hydrogen stable isotope measurements of water by isotope ratio infrared spectroscopy (iris). *Limnology and Oceanography: Methods*, 10(12):1024–1036, 2012. doi: <https://doi.org/10.4319/lom.2012.10.1024>. URL <https://aslopubs.onlinelibrary.wiley.com/doi/abs/10.4319/lom.2012.10.1024>.

Manfred Gröning. Improved water  $\delta^2\text{H}$  and  $\delta^{18}\text{O}$  calibration and calculation of measurement uncertainty using a simple software tool. *Rapid Communications in Mass Spectrometry*, 25(19):2711–2720, 2011. doi: <https://doi.org/10.1002/rcm.5074>. URL <https://analyticalsciencejournals.onlinelibrary.wiley.com/doi/abs/10.1002/rcm.5074>.

Scipy Community. Scipy minimize tutorials. <https://docs.scipy.org/doc/scipy/tutorial/optimize.html#constrained-minimization-of-multivariate-scalar-functions-minimize>. Accessed on 15/12/2022.

Scipy Community. Scipy minimize documentation. <https://docs.scipy.org/doc/scipy/reference/generated/scipy.optimize.minimize.html?highlight=minimize>, b. Accessed on 15/12/2022.

Scipy Community. Scipy linear regression documentation. <https://docs.scipy.org/doc/scipy-1.9.3/reference/generated/scipy.stats.linregress.html#scipy.stats.linregress>, c. Accessed on 15/12/2022.

Wikipedia Contributors. p-value. <https://en.wikipedia.org/wiki/P-value>. Accessed on 31/01/2023.

Boaz Luz and Eugeni Barkan. The isotopic ratios  $^{17}\text{O}/^{16}\text{O}$  and  $^{18}\text{O}/^{16}\text{O}$  in molecular oxygen and their significance in biogeochemistry. *Geochimica et Cosmochimica Acta*, 69(5):1099–1110, 2005. ISSN 0016-7037. doi: <https://doi.org/10.1016/j.gca.2004.09.001>. URL <https://www.sciencedirect.com/science/article/pii/S001670370400660X>.